


ADC Medi-Alert Services, LLC
medi-alert.com

PROGRAMMING A TRANSMITTER INTO THE ADC CONSOLE (ADC 1- 3 AND 4 Models)

PLEASE READ THESE IMPORTANT INSTRUCTIONS

Transmitters such as pendants, wristbands, smoke detectors or motion detectors can be added by following these steps.

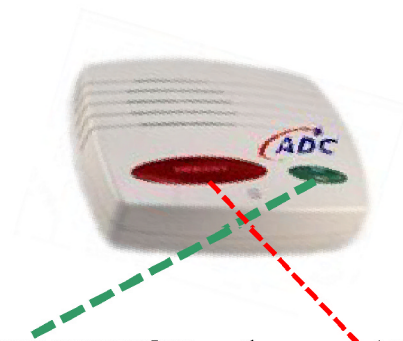

- 1- Press and hold the **GREEN RESET BUTTON** on the ADC Console, after the initial “**BEEP**” upon pressing the **RESET BUTTON** you will hear a second beep “approximately 5 seconds later”. Release the **RESET BUTTON** after the second beep.
- 2- Press the transmitter pendant/wristband and the **RED HELP/EMERGENCY BUTTON** on the ADC Console at the same time. Upon activations, the console will “beep” to let you know the transmitter is loaded.
- 3- Press the **GREEN RESET BUTTON** once and release. The ADC Console is now back in service.